
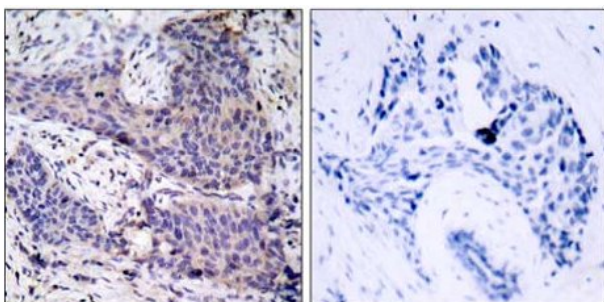


BAD (Ab-136) Antibody

Catalog No.	75384
Product Name	BAD (Ab-136) Antibody
Size	50ug, 100ug
Concentration	1 mg/ml
Specificity	Detects endogenous levels of total BAD protein.
Reactivity	Human, Mouse, Rat
Immunogen	Synthesized non-phosphopeptide derived from mouse BAD around the phosphorylation site of serine136 (S-R-S ^P -A-P).
Applications	WB: 1/500 – 1/2000; IHC: 1/50 – 1/100
Swiss-Prot No.	Q92934 (Human), Q61337 (Mouse), O35147 (Rat)
Host Species	Rabbit
Clonality	Polyclonal
Isotype	IgG
Storage	-20°C for one year. Avoid freeze / thaw cycles.
Purification	Affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.

Images

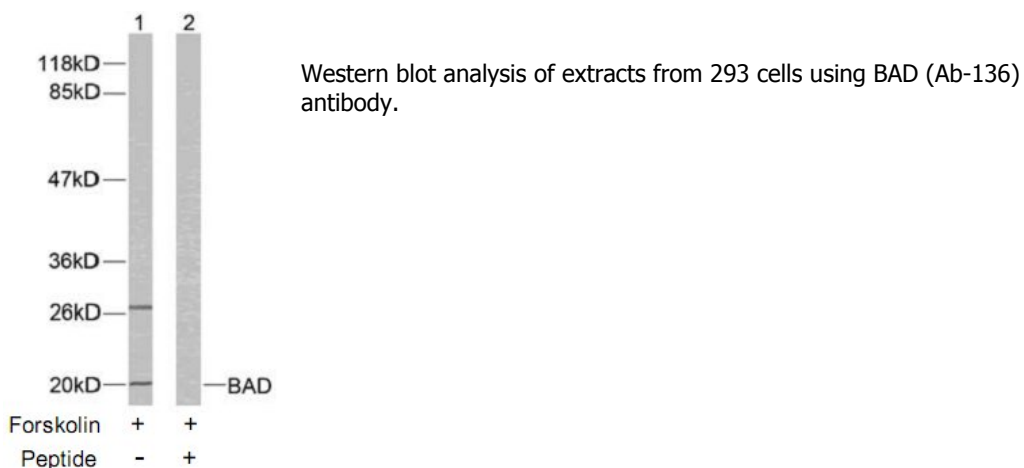

Immunohistochemical analysis of paraffin-embedded human breast carcinoma tissue using BAD (Ab-136) antibody.

Peptide - +

For Research Use Only