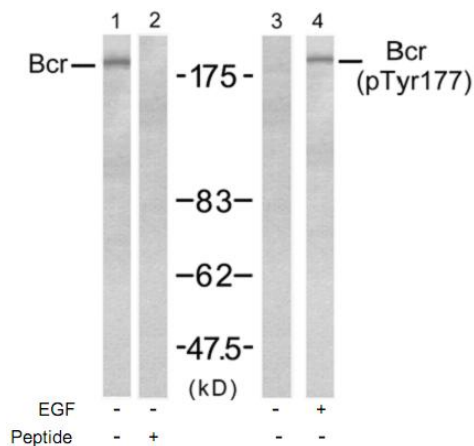


Bcr (Ab-177) Antibody

Catalog No.	75518
Product Name	Bcr (Ab-177) Antibody
Size	50ug, 100ug
Concentration	1 mg/ml
Specificity	Detects endogenous levels of total Bcr protein.
Reactivity	Human, Mouse
Immunogen	Synthesized non-phosphopeptide derived from human Bcr around the phosphorylation site of tyrosine 177 (P-F-Y ^P -V-N).
Applications	Western Blot: 1/500 – 1/1000
Entrez Gene No.	613 (Human); 110279 (Mouse)
Swiss-Prot No.	P11274 (Human)
Host Species	Rabbit
Clonality	Polyclonal
Isotype	IgG
Buffer	Phosphate buffered saline (without Mg ²⁺ and Ca ²⁺), pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol.
Storage	-20°C for one year. Avoid freeze / thaw cycles.
Purification	Affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.

Images

Western blot analysis of extract from A431 cells untreated and treated with EGF (200ng/ml, 5min), using Bcr (Ab-177) antibody (Lane 1 and 2) and Bcr (phospho-Tyr177) antibody (Lane 3 and 4).

For Research Use Only