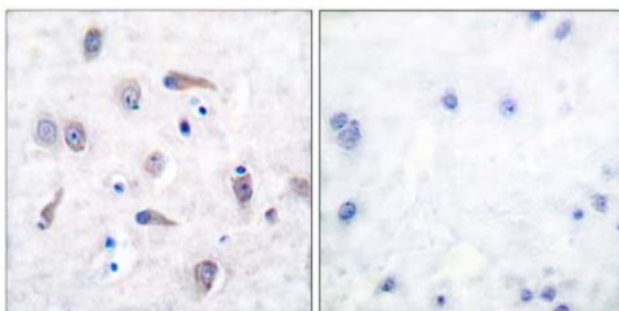



Amyloid β A4 (Ab-743/668) Antibody

Catalog No.	75625
Product Name	Amyloid β A4 (Ab-743/668) Antibody
Size	50ug, 100ug
Concentration	1 mg/ml
Specificity	Detects endogenous levels of total Amyloid β A4 protein.
Reactivity	Human, Mouse, Rat
Immunogen	Synthesized non-phosphopeptide derived from human Amyloid β A4 around the phosphorylation site of threonine743/668 (A-V-T ^P -P-E).
Applications	IHC: 1/50 – 1/100; IF: 1/100 – 1/200
Swiss-Prot No.	P05067 (Human)
Host Species	Rabbit
Clonality	Polyclonal
Isotype	IgG
Buffer	Phosphate buffered saline (without Mg ²⁺ and Ca ²⁺), pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol.
Storage	-20°C for one year. Avoid freeze/thaw cycles.
Purification	Affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.

Images

Immunohistochemical analysis of paraffin-embedded human brain tissue using Amyloid β A4 (Ab-743/668) Antibody.

Left: peptide –
Right: peptide +

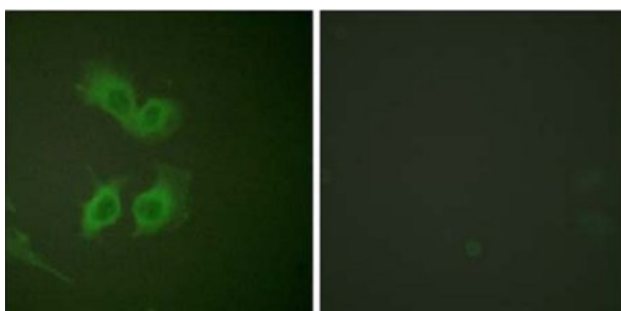

Immunofluorescence analysis of HeLa cells, using Amyloid β A4 (Ab-743/668) antibody.

Left: peptide –
Right: peptide +

For Research Use Only

Full Moon BioSystems, Inc.

754 N. Pastoria Avenue • Sunnyvale, CA 94085

Phone: 408.737.1702 • Toll Free: 877.316.7308 • Fax: 408.732.7230 • www.fullmoonbiosystems.com